

Instrument Processing




Topic: 10



Processing of Items for Reuse

Instruments, Gloves and Glass Syringes

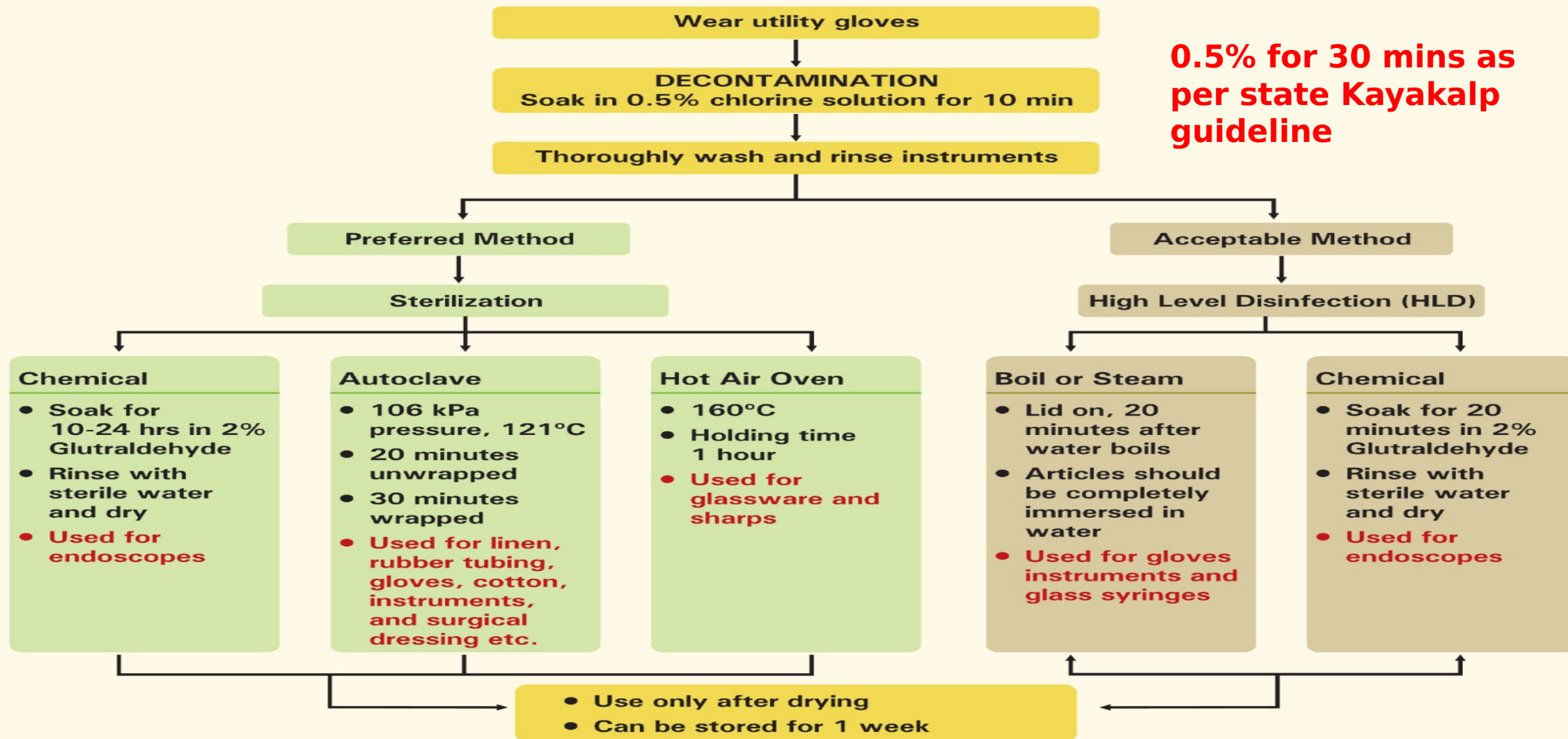

Preparation of 1 Litre Bleaching Solution




Wear utility gloves and plastic apron



- Take 1 L water in a plastic bucket
- Make thick paste in plastic mug with 3 level teaspoons (15 g) bleaching powder and some water from bucket
- Mix paste in water to make 0.5% of chlorine solution

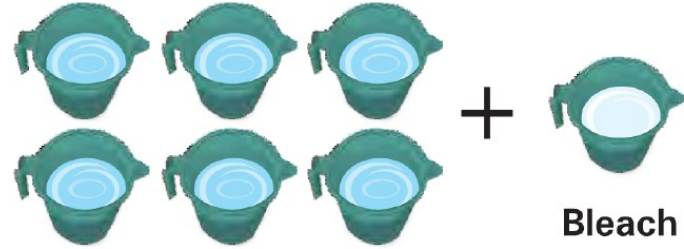

- 3.5%**
- Mix 6 part water with one part of Sodium Hypochlorite solution (Liquid bleach)
- 
- Water + Bleach

- Maintain same ratio for large volumes
- Make fresh solution in every shift and preferably keep covered

Hypochlorite solution available concentration and dilution

Available concentration	Desired Concentration	Dilution
3.5%	0.5%	6 part of water and 1 part of Hypochlorite solution
3.5%	1%	2.5 part of water and 1 part of Hypochlorite solution
5%	0.5%	9 part of water and 1 part of Hypochlorite solution
5%	1%	4 part of water and 1 part of Hypochlorite solution
Bleaching powder	0.5%	1 litre water 15 gms Bleaching powder
Bleaching powder	1 %	1 litre water 30 gms Bleaching powder

High Level disinfection/ Chemical sterilization



BIOMEDICAL WASTE MANAGEMENT RULES, 2016 (AMENDED 2018 & 2019)

Yellow Category

Non-chlorinated plastic bag/container



- Human and animal anatomical wastes
Tissues, organs, body parts and fetus below the viability period
- Soiled waste
Contaminated with blood, body fluids like dressings, swabs, plaster casts, linen, pads, mask and gown
- Blood bags
Pre-treated before putting it in yellow containers
- Chemical wastes
Including discarded disinfectants



- Cytotoxic drugs
- Items contaminated with cytotoxic drugs

Red Category

Non-chlorinated plastic bag/container



- Recyclable waste (plastic)
- Tubing
- Bottles
- IV tube/sets
- Catheters
- Urobags
- Syringes
- Vacutainers
- Gloves

Blue Category

Puncture proof and leak proof container with blue coloured marking



- Waste glassware
- Broken glass
- Medicine vials/ampules
- Metallic body implants

White (Translucent)

Puncture proof, leak proof, tamper proof container



- Waste metal sharps
- Used, discarded and contaminated needles
- Syringes with fixed needles
- Needles cut in hub-cutter scalpel blades